

Product Datasheet

GAPDH Antibody

Catalog No: SB-AB0037 Reactivity: H, M, R, Mk, Ch, Ze, Fish

Isotype: Rabbit IgG

Applications: WB IHC ICC/IF IP FC

ShareBio
圣尔生物

www.share-bio.com

Information

UniProt ID: P04406

All Names: aging-associated gene 9 protein; G3P; G3PD; GAPD; GAPDH; glyceraldehyde 3-phosphate dehydrogenase; Glyceraldehyde-3-phosphate dehydrogenase; MGC88685;

Form: Liquid

Storage instructions: Store at +4°C short term. Store at -20°C long term. Avoid freeze / thaw cycle.

Storage buffer: pH 7.4, 150mM NaCl, 0.02% sodium azide and 50% glycerol.

Purity: Affinity-chromatography

Immunogen: A synthesized peptide

Molecular Wt.: 37 kDa

Application

WB: 1:3000~1:10000

IHC: 1:100~1:500


IF/ICC: 1:100~1:250

IP: 1:50

FC: 1:50

Background

Glyceraldehyde 3 phosphate dehydrogenase (GAPDH) is well known as one of the key enzymes involved in glycolysis. GAPDH is constitutively abundant expressed in almost cell types at high levels, therefore antibodies against GAPDH are useful as loading controls for Western Blotting. Some pathology factors, such as hypoxia and diabetes, increased or decreased GAPDH expression in certain cell types.

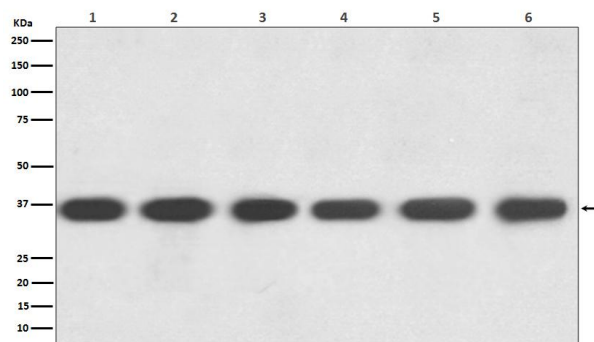

Western blot analysis of (1) Hela, (2) Jurkat whole cell lysates, (3) Mouse kidney lysates, (4) Mouse spleen lysates, (5) RAW 264.7 whole cell lysates, (6) Rat brain lysates with GAPDH Antibody.

For Research Use Only. Not For Use In Diagnostic Procedures.

www.share-bio.com